




Sterimed Oxygen Mask

SMD 704

-  Oxygen Mask
-  Máscara De Oxígeno
-  Máscara de oxigênio
-  Masque d'oxygène

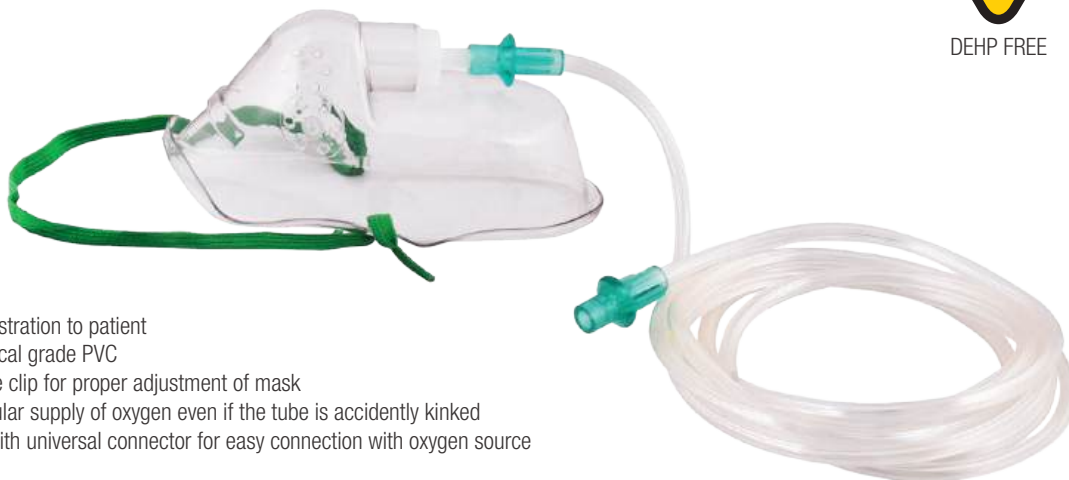

DEHP FREE

- Used for efficient oxygen administration to patient
- Made from soft, non-toxic, medical grade PVC
- Elastic strap and aluminum nose clip for proper adjustment of mask
- Multi-channel tube ensures regular supply of oxygen even if the tube is accidentally kinked
- Proximal end of tube provided with universal connector for easy connection with oxygen source

Packing

- Non- Sterile, individually packed in a ribbon.
- Box of 20; Master Carton of 200.



Gross: 18.4kg (Adult)



740 x 425 x 425mm

SIZES	TUBE LENGTH
Adult	210cm (7ft)
Child	210cm (7ft)



Sterimed Non-Rebreather Mask

High Concentration Non-Rebreather Oxygen Mask

SMD 704 N

-  Masker Non-Rebreather Sterimed
-  Máscara sin rebreather Sterimed
-  Máscara Sterimed sem reinalação
-  Masque sans réinspiration Sterimed

- Specially designed for effective high-concentration oxygen therapy.
- Made from medical-grade, clear, odorless PVC for patient comfort and visibility.
- Atraumatic feathered edges reduce skin pressure and breakdown risk.
- Contoured mask seals naso-oral area, reducing ambient air entry.
- Adjustable elastic strap for secure, customizable fit.
- Reservoir bag stores 100% oxygen for consistent high FiO₂.
- Design supports variable tidal volumes and breathing rates.
- Unidirectional valves prevent rebreathing and CO₂ retention.
- Lightweight reservoir improves comfort and oxygen flow.

Packing

- Non-Sterile, individually packed.
- Box of 1, Master Carton of 100.



Gross: 10.0kg



670 x 300 x 480mm

SIZES	TUBE LENGTH
Adult	200cm
Child	200cm