






Catheter Mount

SS 758

-  Kateter Mount
-  Montaje del catéter
-  Montagem de cateter
-  Support de cathéter

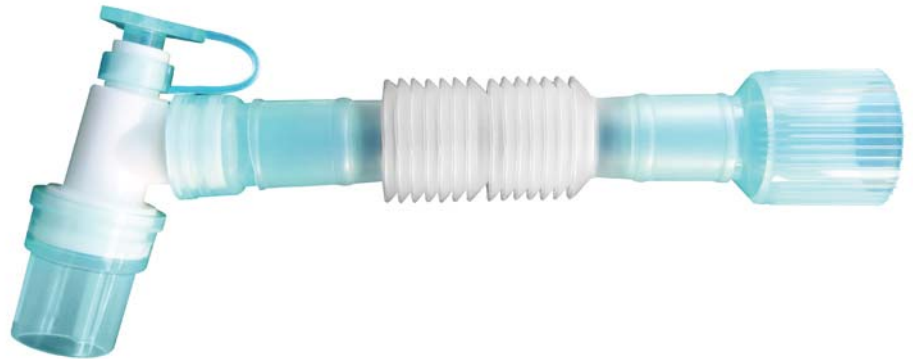

- Provided with Suction & Bronchoscopy Ports.
- Compatible with all Breathing and Ventilator circuits.
- Double swivel connector and flexible tubing provides mobility and flexibility to the patient end of the circuit.
- The tubes are made of a special flexible co-polymer and have a minimal compliance with low dead space volume.
- The catheter mount is fitted with a double swivel connector
- Each catheter mount has a 22 mm female / 15 mm male end to connect to the breathing system and 15 mm female end to connect to the endotracheal tube or laryngeal mask.

Packing

- EO Sterile/ Disposable/ Individual Packed.
- Master Carton of 100 pcs.



Gross: 5.5kg



290 x 280 x 330mm

HME Filter

SS 752

-  HME Filter
-  HME Filtro
-  Filtro HME
-  HME Filtre



Packing

- Individually Blister Packed.
- Master Carton of 100 pcs.



Gross: 6kg



430 x 380 x 330mm

- Used in mechanically ventilated and laryngectomy patients to heat and moisturize the air and help prevent complications, that arise when patients lose the ability to breathe through their nose and upper airway. Captures heat and moisture on expiration and returns it to the patient in inspiration
- Suitable for spontaneously breathing patients with a tracheostomy tube
- Electrostatic filter with 99.99% filtration capacity
- Lightweight, compact design reduces circuit weight
- Low resistance to flow reduces work of breathing
- ISO standard 15 mm and 22 mm fitting connects to breathing circuit
- Supplied clean and individually packaged
- Available in adult and paediatric options