

Reinforced Endotracheal Tube - Cuffed & Plain

(Nasal/Oral)

Cuffed | SMD 719 C

Plain | SMD 719 P

Reinforced Endotracheal Tube- Cuffed & Plain


Tubo Endotraqueal Reforzado Plano

Endotraqueal reforçado Tube - Cuffed & Plain

Endotrachéale renforcé Tube - Revers & Plain



PHTHALATE FREE

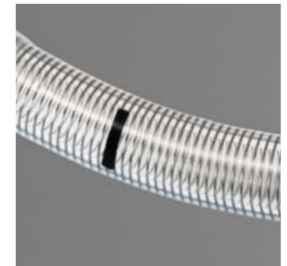

Uniform spiral metal to prevent tube kinking

- Used to maintain a patent airway and mechanically ventilate patients
- Made from non-toxic, non-irritant medical grade PVC
- Uniform spiral metal-reinforced tube wall to prevent the tube from kinking
- Reinforcement supports tube position and verification during X-ray
- Accurate and safe positioning ensured by graduation marks and tip-to-tip radiopaque line
- High volume / low pressure cuff ensures efficient low pressure cuff seal (For cuffed tubes)
- Universal 15 mm connector at proximal end
- Recommended use of bite block / Guedel airway
- E.T.O sterile and pyrogen free. For single patient use

Packing

- Sterile, Individually packed in peelable pouch / blister packed.
- Box of 25; Master Carton of 200 pcs.



Gross: 6.7kg (300 pcs)
Gross: 12.5kg (300 pcs)



660 x 320 x 380mm (3.0 to 4.5mm)
660 x 380 x 500mm (5.0 to 9.5mm)

SIZES - CUFFED

3.0 mm, 3.5 mm, 4.0 mm, 4.5 mm, 5.0 mm, 5.5 mm, 6.0 mm, 6.5 mm, 7.0 mm, 7.5 mm, 8.0 mm, 8.5 mm, 9.0 mm, 9.5 mm

SIZES - PLAIN

3.0 mm, 3.5 mm, 4.0 mm, 4.5 mm, 5.0 mm, 5.5 mm, 6.0 mm, 6.5 mm, 7.0 mm, 7.5 mm, 8.0 mm, 8.5 mm